


DDC Raw Capture

Instructions to generate DDC Raw Capture images:

1. File -> Open Archives -> Select archives to open. You may select more than one.
2. Open DDC Raw viewer if not already open: Data Viewers -> DDC Raw
3. Maximize the window
4. Set settings on the right side:
 - CCD ID: Red or BLUE
 - Chart Height: 500
 - Sync Axes: Both
 - Show Regions: checked
 - Line Window: 1
 - Right click the "All Taps Video" chart's y-axis, set range [14500, 16500]
 - Right click the "All Taps CDS" chart's y-axis, set range [-50, 500]
 - If desired, right click on each chart and select "All series data points visible"
 - All other settings can be left at the defaults
5. Capture screenshot of whole window, except stop before the "All Taps Raw Samples" chart. Include the settings area on the right side.
6. Repeat as needed for RED/BLUE

Below is an example screenshot:

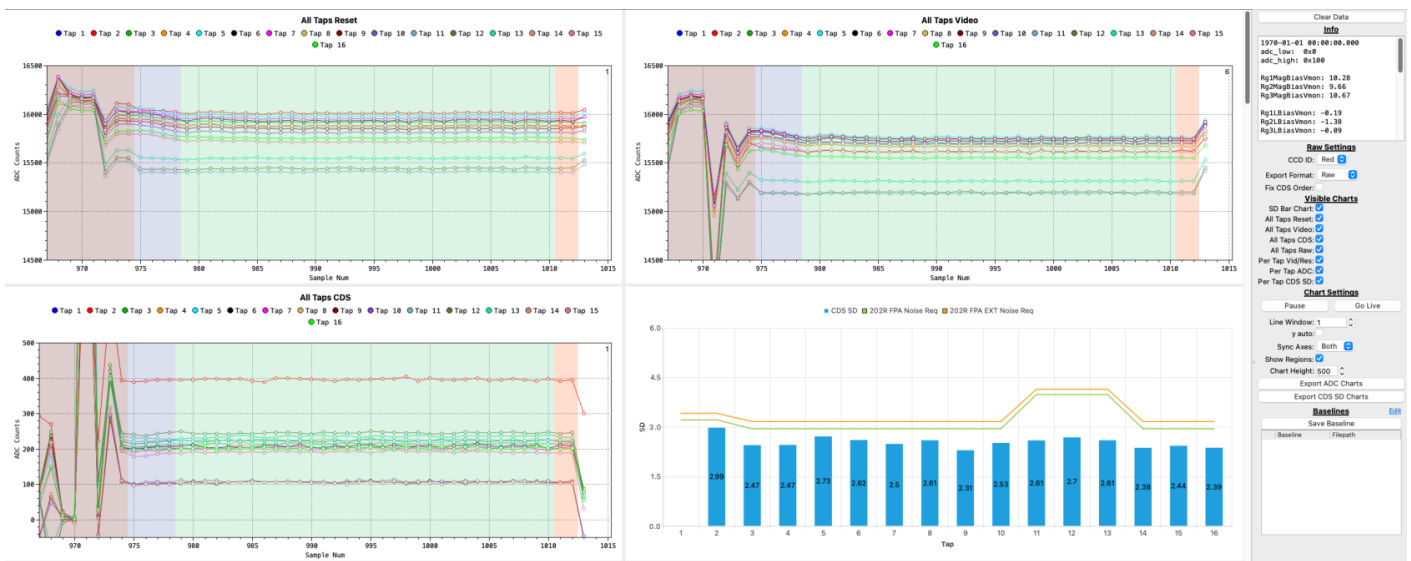

Revision #9

Created 22 February 2024 17:41:48 by Bradley Tse

Updated 23 February 2024 15:26:06 by Bradley Tse